


**SURGEINE**<sup>®</sup>  
*Surely Hygiene*

STERILIZATION WRAP

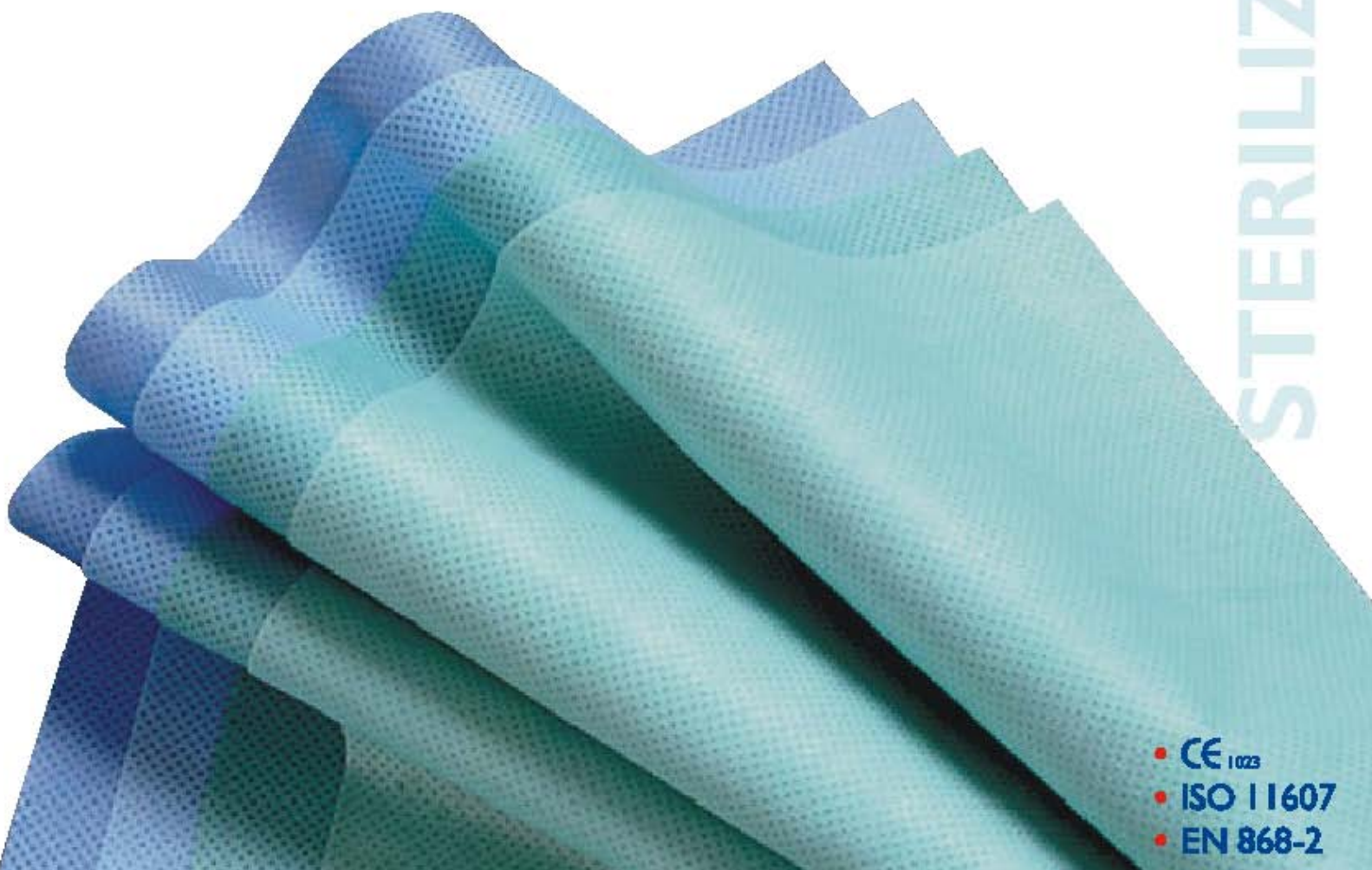

- CE<sub>1023</sub>
- ISO 11607
- EN 868-2

**SURGEINE** sterilization wraps are made from SMMS material consist of spunbond/melt-blown/melt-blown/spunbond. “spunbond” provides the strength and “meltblown” provides the barrier. These multiple layer provides excellent protection from microbial contamination and gives adequate strength.

Surgeine sterilization wraps conforms to international standard ISO 11607-1:2006 also have CE mark and passes all the tests of EN 868-2.

## Characteristics of **SURGEINE** Sterilization Wrap

### **BARRIER EFFECTIVENESS**

Prevent microbial penetration and maintain sterility of surgical pack and prevent penetration of liquids.

### **PENETRABILITY & AERATION**

Allows steam, EO and Plasma to penetrate & aeration post-sterilization

### **DRAPEABILITY**

Easy to use, very soft, free from creases, conforms to equipment pack contours smoothly and closely.

### **STRENGTH**

Puncture and tear resistance, adequate strength than crepe.

### **LINTING**

Naturally low linting to help maintain the integrity of sterile pack.

### **IDENTIFICATION**

Each sheet is provided with individual self adhesive label to write pack's internal contents. This will help you to verify pack's internal contents before opening the pack.

Apart from that surgeine sterilization wrap is anti static, too flexible, non reactive, non-toxic and low cost in use.

Surgeine offers three different types of sterilization wraps. You can choose from our wide range to suits your requirements of wrapping strength and sizes.

Our range consists of :

## **REGULAR WRAPS**

Single sheet of SMMS non woven fabric, available in different colors. You can use a single sheet for regular wrapping or two sheets sequentially wrapped one after another of different color for double wrapping. Our regular strength wraps are ideal for light linen packs, gowns, single basins and another light to medium weight and size of surgical packs and trays.

## **HEAVY DUTY WRAPS**

Our heavy duty wraps are constructed of two sheets of SMMS fabric. These two Sheets are bonded together on there two side edges with the help of ultrasonic. This wrap gives you extra strength for your surgical packs and trays also saves your time. This wraps is suitable for Basin sets, Instrument sets, procedure trays and large linen packs.

## **SUPER HEAVY DUTY WRAPS**

These wraps are constructed of three sheets of non woven fabric. Two layers of SMMS non woven fabric and one layer of Spun bond non woven fabric in the middle are bonded together on there two side edges with the help of ultrasonic. Super heavy duty wraps are ideal for packaging heavy instruments sets, reactors, drills and another very heavy weight and size of surgical packs and trays.

## REGULAR WRAPS

Item #	Size	Quantity
2001	12" x 12"	1000/cs
2002	18" x 18"	1000/cs
2003	24" x 24"	1000/cs
2004	36" x 36"	500/cs
2005	48" x 48"	500/cs
2006	54" x 54"	300/cs

## HEAVY DUTY WRAPS

Item #	Size	Quantity
2011	18" x 18"	1000/cs
2012	24" x 24"	1000/cs
2013	36" x 36"	500/cs
2014	48" x 48"	500/cs
2015	54" x 54"	300/cs
2016	54" x 72"	100/cs

## SUPER HEAVY DUTY WRAPS

Item #	Size	Quantity
2021	24" x 24"	1000/cs
2022	36" x 36"	500/cs
2023	48" x 48"	500/cs
2024	54" x 54"	300/cs
2025	54" x 72"	100/cs
2026	54" x 90"	100/cs

**SURGEINE HEALTHCARE (INDIA) PVT. LTD.**  
 101, Vardhman Tri Plaza, CU Block, Pitampura, New Delhi-110088 (India)  
 Tel. : +91-11-47552462, 27348301 Telefax : +91-11-47552424  
 Mob. : +91-9350990935, 91-9310990931  
 E-mail : sales@surgeine.com

- CE<sub>1023</sub>
- ISO 11607
- EN 868-2
- ISO 9001
- ISO 13485
- GMP